



Medical Technologies


EasyOne Pro/LAB

Complete point-of-care pulmonary function testing

Making space for what matters in pulmonary health

The **EasyOne Pro/LAB** revolutionizes pulmonary function testing with its comprehensive range of tests, compact design, and user-friendly interface. Its capabilities, ranging from spirometry to lung volumes and DLCO, cover a wide spectrum of diagnostic needs in respiratory health.

The **EasyOne Pro/LAB** enables healthcare professionals to efficiently evaluate lung function, aiding in diagnosis, monitoring, and management of respiratory conditions. The combination of portability, ease of use, and accuracy makes the **EasyOne Pro/LAB** suitable for various clinical settings.

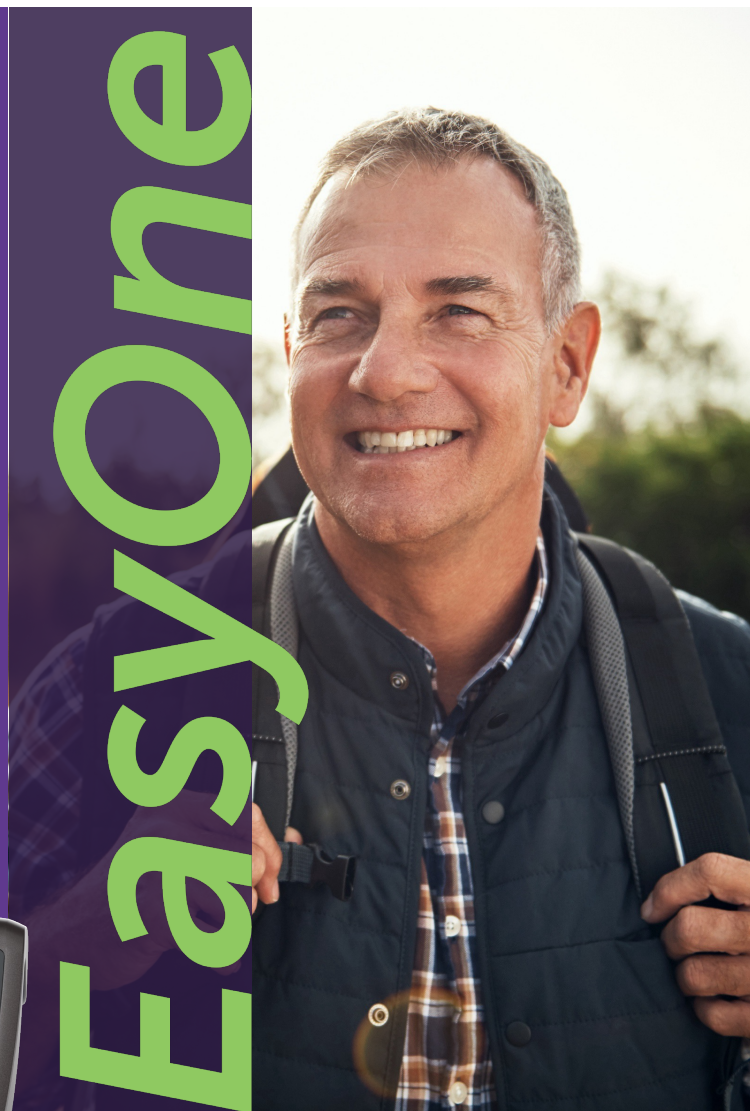

Precise performance



EasyOne Pro

- Spirometry
 - DLCO
 - Provocation testing
- Lung volumes**
- Single-breath helium dilution

EasyOne Pro LAB

- Spirometry
 - DLCO
 - Provocation testing
- Lung volumes**
- Multiple-breath nitrogen washout
 - Single-breath helium dilution



Specifications

- EasyOne Pro:** FVC, FVL, SVC, MVV, Provocation and Pre/Post testing, DLCO
- EasyOne Pro LAB:** FVC, FVL, SVC, MVV, Provocation and Pre/Post testing, DLCO, MBW, LCI
- Languages:** Chinese, Croatian, Danish, Dutch, English, Finnish, French, German, Italian, Japanese, Norwegian, Polish, Portuguese, Russian, Spanish, Swedish, Turkish, Vietnamese
- EMR connectivity:** HL7, XML, GDT
- Dimensions:** 10.6 x 13.2 x 10.6", 17.65 lb
27 x 33.5 x 27 cm, 8 kg
- Spirometry Standards:** ATS/ERS 2022, 2019 & 2005, NIOSH, OSHA, SSA Disability
- DLCO Standards:** ATS/ERS 2017 & 2005

Article numbers

- 3000-1:** EasyOne Pro
- 3100-1:** EasyOne Pro LAB
- 2050-1:** Spirette box of 50
- 2050-5:** Spirette box of 200
- 2050-10:** Spirette box of 500
Not available in all countries
- 3050-1:** DLCO Barriettes box of 50
- 3050-2:** DLCO Barriettes box of 100
- 3150-1:** FRC Barriettes box of 40 (for use with EasyOne Pro LAB only)
- 3150-2:** FRC Barriettes box of 80 (for use with EasyOne Pro LAB only)



ndd Medical Technologies
300 Brickstone Square, Suite 604 Andover, MA, 01810 USA
Phone: +1-978-470-0923

ndd Medizintechnik AG
Technoparkstrasse 1, CH-8005 Zurich, Switzerland
Phone: +41 44 512 65 00



Features

- Portable and lightweight design tailored for testing at the point of care.
- User-centric approach for intuitive operation, instant usability with no warm-up time.
- Hygienic solution with Spirette and Barriette consumables eliminates the risk of cross-contamination.
- Swiss precision engineering
- Minimal maintenance, allowing more time to be dedicated to patient care.

TrueFlow

makes the difference

ndd's unique ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity.

ndd **TrueFlow** is a resistance-free solution that does not require calibration during its lifetime.

TrueCheck

automated precision

- Patented method to verify performance of ndd machines
- Eliminates the need for gas-calibration as well as the need for an external DLCO simulator
 - 5-point calibration including stability and linearity of tracer gas and CO sensors

EasyOne Connect

EasyOne Connect software offers one database and one platform for EasyOne point-of-care-solutions. Data management has never been easier.

