




Medical Technologies

Get to know Buddy,
your spirometry coach.

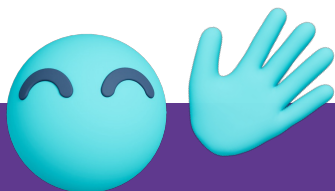

The future of spirometry

The **EasyOne Sky** transforms respiratory diagnosis and management with a fully connected spirometry solution.

Guided by Buddy, n d d's built-in spirometry coach, patients can perform tests more easily while providers achieve consistent, high-quality results. Our ecosystem supports healthcare systems at every stage, from deployment through diagnosis. Secure cloud infrastructure protects every test, while real-time insights provide visibility into testing performance helping healthcare teams deliver better care with confidence.



EasyOne Sky

Guided spirometry. Connected insights.


EasyOne



Spirometry designed for your workflow.



Hi, I am Buddy,
your spirometry
coach!

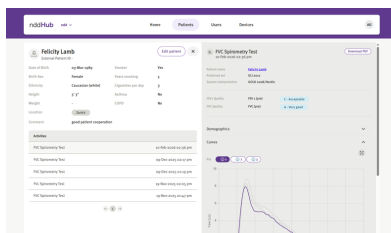

- Buddy, your built-in spirometry coach, guides patients through each maneuver and provides real-time feedback to help achieve high-quality test results.
- Powered by nddCloud, **nddHub** is a web-based portal that provides centralized fleet and performance management across your entire network.
- Standardized setup, EMR integration, and structured support enable smooth spirometry rollout across any care setting.
- RT-led training and continuous education keep clinical teams confident and testing consistent over time.
- HITRUST-certified infrastructure automatically backs up and protects every test while enabling secure access anywhere.

Specifications

Test types: FVC, FVL, Pre/Post-testing

Languages: English, Spanish

Dimensions: 4.3" x 3.3" x 1.3", 4.9 oz / 111 x 84 x 33 mm, 140 g

Parameters: BEV, FEF₂₅, FEF₂₅₋₇₅, FEF₅₀, FEF₇₅, FET, FEV_{0.75}, FEV_{0.75}/FVC, FEV₁, FEV₁/FVC, FEV₁Q, FIVC, FVC, PEF, PIF, FVC %change, FEV1 %change

Predicted normal values spirometry: GLI (Quanjer) 2012, NHANES-III (Hankinson) 1999, Global GLI (Bowerman) 2023

Standards: ATS/ERS 2019 Spirometry Standard, GOLD 2026

Data management: EasyOne Mobile nddCloud nddHub

App compatibility: Windows Download on the App Store GET IT ON Google Play
iOS and Android not available in all countries

Article numbers

9500-1 EasyOne Sky Comfort package

9500-2 EasyOne Sky Enterprise package

5050-50 EasyOne FlowTube box of 50

5050-200 EasyOne FlowTube box of 200



ndd Medical Technologies
300 Brickstone Square, Suite 604 Andover, MA, 01810 USA
Phone: +1 978 470 0923

ndd Medizintechnik AG
Technoparkstrasse 1, CH-8005 Zurich, Switzerland
Phone: +41 44 512 65 00



Advancing respiratory care through connected data.

- Real-time PFT testing insights
- Comprehensive fleet management
- User & clinic performance visibility
- Centralized healthcare network oversight



ndd's unique ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity.

ndd **TrueFlow** is a hygienic and resistance-free solution that does not require calibration during its lifetime.



ndd's modern and secure connectivity solution is designed to connect to healthcare platforms and EHR/EMR systems. It is built on state-of-the-art cybersecurity standards.

